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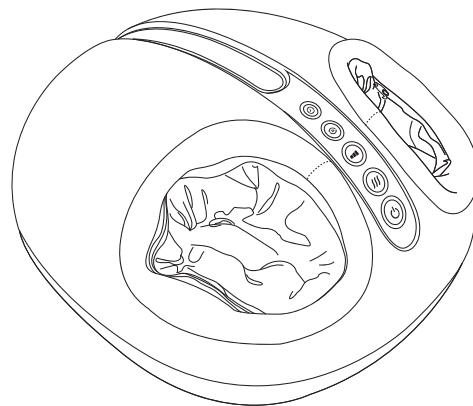
Foot Massager

USER MANUAL

Thank you for purchasing and
using our products

Read through the entire manual
before installation and use!


FOOT MASSAGER SHIATSU PREMIUM

Read the user manual carefully to make sure that the product is assembled, used and maintained correctly, as described in the instructions in this user manual. Keep the user manual for future reference.

SAFETY PRECAUTIONS

- In order to ensure safe & proper use of the product, please refer to the manual before use.
- This product can be used by children of 8 years of age and older as well as by people with reduced physical, sensory and mental capacity or lack of experience and knowledge under supervision or if they have been given instructions relating to use of this product in a safe manner and understand potential risks.
- Supervise children, to ensure that they do not play with this product.
- This product has a heater which means that it can become very hot. If you are sensitive to heat you must pay attention during use to avoid discomfort.
- If the power cord is damaged it must be replaced by the manufacturer or other similarly qualified person.
- Do not try to repair the cord or any other parts by yourself.
- Keep this product away from sharp objects and heat sources.
- If the cloth cover is broken, stop using the product.
- After use, always pull out the power plug from the socket.
- Before cleaning and maintenance, always pull out the power plug from the socket.
- Do not cover the unit while it is on. This can cause failure, damage or over heating on the motor unit.
- If you feel any discomfort during use, immediately stop using the product.
- If the product malfunctions or is damaged, please do not operate. Contact the manufacturer for inspection and reparation.
- Do not use accessories not provided by the manufacturer.
- Do not operate the product with wet hands.
- Do not let water or other liquids touch the controls.
- Do not use the unit in a very humid or dusty environment.
- Do not use the unit in places where the temperature exceeds +40 °C.
- For indoor use only.
- Always consult your doctor before use if:
 - You are under treatment.
 - You have high blood pressure, heart disease, cerebral arteriosclerosis, cerebral thrombosis, malignant tumor, bone osteoporosis.
 - You are pregnant.
 - You are using a pacemaker or other implanted medical apparatus.
- Do not use this appliance while you are asleep.
- Do not use if you are not feeling well.
- Do not use for longer than 30 minutes at a time.

DESCRIPTION

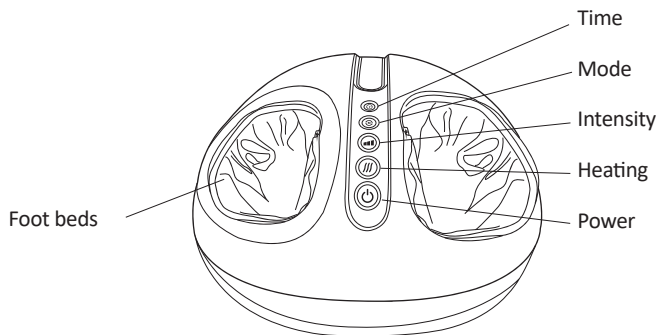WARNING!

If these instructions are not followed, an accident may lead to serious personal injury or property damage.

OPERATION MODE

TIME KEY: To control the massage time

Blue color: Set time is 10 minutes.

Red color: Set time is 15 minutes.

Purple color: Set time is 20 minutes.

I as the defaults mode, first press enter into mode II, second press change into mode II.

Mode I: blue lights, mode II: red light, mode II1: flashing alternating red and blue lights.

Mode I: soft sole massage + pneumatic massage.

Mode II: pneumatic massage only.

Mode III: comfortable sole massage + pneumatic massage.

INTENSITY KEY: Air pressure intensity control

For low intensity: Blue light

For medium intensity: Green light

For high intensity: Red light

HEATING KEY: Heating control function

Press the key to turn off or turn on the heat

POWER KEY: turn off or turn on

NOTE! This foot massager is equipped with an overheating protection. If the temperature exceeds the safety temperature, it will automatically stop working. If this happens, shut off the power and disconnect the plug from the socket. Allow the appliance to rest 15 minutes before using again.

CLEANING AND MAINTENANCE

STORAGE

When not in use, put the Foot Massager Shiatsu Premium in its box and store it in a dry, cool place. Make sure not to tweak the cord. Avoid sharp object so that the surface is not scraped or broken. Do not store in direct sunlight or in a high temperature. If you will not be using the product for a long time, use a dust cover to prevent dust from getting in to the massager.

CLEANING

- Before cleaning, make sure that the power is off and that the plug is not connected. Use a dry cloth to clean the surface of the massager. Do not use a corrosive cleaner.
- The fabric cover can be removed and cleaned separately by hand wash only.

TROUBLESHOOTING

If you experience any problems during use, please refer to the following tips. If the product still doesn't work, please contact customer service at Rusta.

| PROBLEM | CAUSE | SOLUTION |
|-------------------------|---|--|
| Massager does not work. | No power. | Access the massager to a power supply.. |
| A sudden stop. | Power is off. | Check that the power supply is connected.. |
| | Massage time has exceeded 20 minutes. | Re-start by pressing the on/off button. |
| | Working time is too long, the massager is overheated. | Re-start after unplugging the massager and letting it rest for at least 15 minutes. |
| | Massage heads are blocked. | Make sure that your feet are relaxed during massage. Make sure that nothing else is trapped in the massager. |

TECHNICAL DATA

| | |
|----------------|---------------|
| Rated voltage | 12V === |
| Rated duration | 10-20 minutes |
| Rated power | 36W |

RIGHT OF COMPLAINT

By law the product must be returned to the place of purchase along with the original receipt in the event of a complaint. The user is liable for any damage to the product that is caused by not using the product for its intended purpose or by not following this user manual correctly. The right of complaint will not apply in these instances.

INFORMATION ABOUT HAZARDOUS WASTE
Used products should be handled separately from household waste. They should be taken for recycling in accordance with local rules for waste management.